



FOCUS ON PRECISION

FOCUS ON INNOVATION

FOCUS ON QUALITY



FOCUS ON
SURGICAL LIQUIDS

EFTIAR® OCTANE/DECALIN

EFTIAR® OFFERS THE VITREORETINAL SURGEON A TEMPORARY TAMPONADE FOR EFFICIENT REATTACHMENT OF THE RETINA, RELEASE OF SUBRETINAL FLUID, DISSECTION OF FIBRO-VASCULAR PROLIFERATION AND CONTROLLING RETINAL BLEEDING.

EFTIAR® HAS CLEAR ADVANTAGES:

- NATURALLY PURE
- 100% FLUORINATED AND PURIFIED
- CHEMICALLY INERT LIQUID
- TRANSPARENT, COLORLESS AND OPTICALLY CLEAR
- LOW REFRACTIVE INDEX
- CLEAR DIFFERENTIATION HIGH SPECIFIC GRAVITY
- SMOOTH FLATTENING OF A DETACHED RETINA
- PASSIVE DRAINAGE OF SUBRETINAL FLUID

EFT.05OCT



EFTIAR® PERFLUORO-OCTANE, 5 ML. VIAL (BOX/1, STERILE)


EFT.0.7OCT



EFTIAR PERFLUORO-OCTANE, 7 ML., VIAL (BOX/1, STERILE)

EFTIAR® OCTANE

- HIGH VAPOR PRESSURE ALLOWS EVAPORATION OF RESIDUAL DROPLETS
- SPECIFIC GRAVITY: 1,77 G/CM3 AT 20°C
- REFRACTIVE INDEX: 1,27 AT 20°C
- VISCOSITY: 1,4 MPAS AT 25°C
- BOILING POINT: 101,3-106°C
- PACKAGING: 1 VIAL OF 5 ML

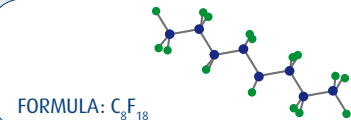

EFT.05DEC



EFTIAR® PERFLUORO-DECALINE, 5 ML. VIAL (BOX/1, STERILE)


EFT.0.7DEC



EFTIAR PERFLUORO-DECALIN, 7 ML., VIAL (BOX/1, STERILE)

EFTIAR® DECALIN



- HIGH INTERFACIAL TENSION TO SILICONE OIL
- SPECIFIC GRAVITY: 1,93 G/CM3 AT 20°C
- REFRACTIVE INDEX: 1,31 AT 20°C
- VISCOSITY: 5,53 MPAS AT 25°C
- BOILING POINT: 140,4-142,4°C
- PACKAGING: 1 VIAL OF 5 ML
- EFT-05-OCT: 5 ML VIAL EFT-05-DEC: 5 ML VIAL

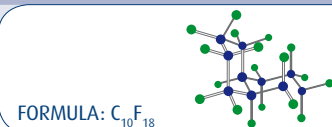

INDICATIONS FOR USE

EFTIAR OCTANE® AND EFTIAR DECALIN® LIQUID IS INDICATED FOR USE AS AN AID IN TEMPORARY FLATTENING OF THE RETINA. IT IS ALSO INDICATED AS AN AID IN MANIPULATION OF THE RETINA DURING SURGICAL TREATMENT OF A RETINAL DETACHMENT. IN PARTICULAR, EFTIAR MAY BE USED IN THE SURGICAL MANAGEMENT OF:

- RETINAL DETACHMENT WITH GIANT TEARS
- RETINAL DETACHMENT WITH PROLIFERATIVE VITREORETINOPATHY
- TRAUMATIC RETINAL DETACHMENT
- SUBLUXATED LENSES
- RETINAL BLEEDING

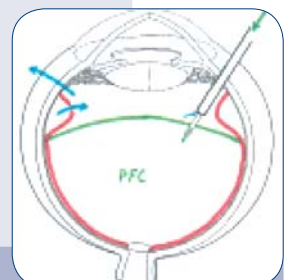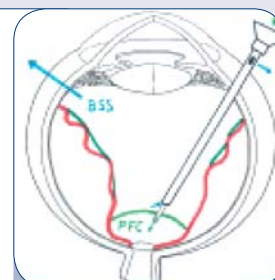

RETINAL ATTACHMENT WITH EFTIAR THROUGH EFTIAR DUAL BORE NEEDLE

SIL-1000®/SIL-5000®

THE D.O.R.C. SILICONE OIL IS AN ULTRA PURIFIED SILICONE OIL WHICH LEADS A MAXIMUM INTERFACIAL TENSION AND MINIMIZES INTERACTIONS BETWEEN TISSUES, CELLS AND ENDO-TAMPONADES MEDIA. THE PHYSICAL PROPERTIES INCLUDE A COMBINATION OF SPECIFIC GRAVITY, REFRACTIVE INDEX AND SURFACE TENSION. THE CHOICE OF VISCOSITY OFFERS AN OPTIMUM BALANCE BETWEEN EASY INJECTION AND A STABLE TEMPORARY TAMPONADE.

SIL.1000



SILICONE OIL, 1,000 CSTS, 10 ML. VIAL (BOX/1, STERILE)


SIL.1000.S



SILICONE OIL, 1,000 CSTS, 10 ML. SYRINGE (BOX/1, STERILE)

SIL-1000® SILICONE OIL

- VISCOSITY: 1.000-1500MPAS
- REFRACTIVE INDEX: 1,40
- SPECIFIC GRAVITY: 0,97 G/CM3 AT 25°C
- SURFACE TENSION: 21 MN/M AGAINST AIR
- INTERFACIAL TENSION: 40 MN/M AGAINST WATER
- PACKAGING: - SIL-1000: 1 VIAL OF 10 ML
- SIL-1000-S: 1 SYRINGE OF 10 ML

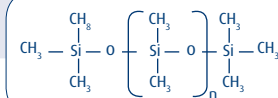

SIL.5000



SILICONE OIL, 5,000 CSTS, 10 ML. VIAL (BOX/1, STERILE)


SIL.5000.S



SILICONE OIL, 5,000 CSTS, 10 ML. SYRINGE (BOX/1, STERILE)

SIL-5000® SILICONE OIL

- VISCOSITY: 5.000-5900MPAS
- REFRACTIVE INDEX: 1,40
- SPECIFIC GRAVITY: 0,97 G/CM3 AT 25°C
- SURFACE TENSION: 21 MN/M AGAINST AIR
- INTERFACIAL TENSION: 40 MN/M AGAINST WATER
- PACKAGING: - SIL-5000 : 1 VIAL OF 10 ML
- SIL-5000-S: 1 SYRINGE OF 10 ML


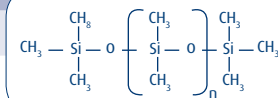


INDICATIONS FOR USE

D.O.R.C. SIL-1000/5000 IS INDICATED FOR USE AS A SHORT TERM TAMPONADE FOR RE-ATTACHMENT OF THE RETINA. IN PARTICULAR, THE PRIMARY INDICATIONS FOR SIL-1000/5000 ARE THE SURGICAL MANAGEMENT OF:

- PROLIFERATIVE VITREORETINOPATHY (PVR)
- RETINAL DETACHMENT CAUSED BY TRAUMA
- GIANT TEARS
- DIABETIC TRACTION DETACHMENT
- CYTOMEGALOVIRUS (CMV) RETINITIS

A SECONDARY INDICATION OF SIL-1000® AND SIL-5000® IS THE SURGICAL MANAGEMENT OF ROP RELATED TRACTION DETACHMENTS.

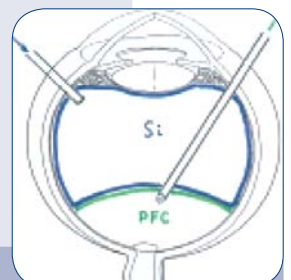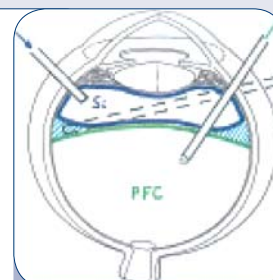

EFTIAR / SIL-1000 EXCHANGE PROCEDURE WITH PASSIVE ASPIRATION THROUGH DORC BACKFLUSH BRUSH NEEDLE

ACCESSORIES FOR OPHTHALMIC LIQUIDS

ACCESSORIES FOR EFTIAR

EFD.100



DISPOSABLE DUAL BORE PFC CANNULA FOR SIMULTANEOUS INFUSION OF HEAVY LIQUIDS AND ASPIRATION OF INTRA-OCULAR FLUIDS. DUAL BORE ENABLES SURGEON TO CONTROL CONSTANT INTRAOCULAR PRESSURE DURING INJECTION, STERILE, BOX/5.

EFT.100



DUAL BORE CANNULA 20 GAUGE CANNULA WITH 25 GAUGE INNER CANNULA FOR SIMULTANEOUS INFUSION OF HEAVY LIQUIDS PERFLUORCARBONS) AND ASPIRATION OF INTRAOCULAR FLUIDS.

EFD.06



EFTIAR® DUAL BORE CANNULA, 23 GAUGE/0.6 mm (BOX/5, STERILE)


ACCESSORIES FOR VFI

1363.DD



VFI PACK FOR D.O.R.C. SILICONE OIL SYRINGES READY TO USE VISCOUS FLUID INJECTION PACK, STERILE, BOX/5.

1363.D

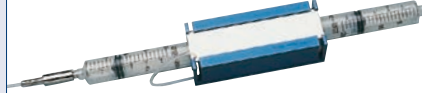

V.F.I. DUAL SYRINGE HOLDER FOR VIALS FOR VISCOUS FLUID INJECTION. HOLDS TWO DISPOSABLE 20CC SYRINGES. COMPLETE WITH PRESSURE TUBINGS, STAINLESS STEEL CONNECTOR AND 40 DISPOSABLE SYRINGES. COMPATIBLE WITH D.O.R.C. SYSTEMS.

1365 D.O.R.C.® VFI ADAPTOR

FOR CONNECTION BETWEEN D.O.R.C.® UNITS AND 1363-DD.

1366 MILLENNIUM® VFI ADAPTOR

FOR CONNECTION BETWEEN MILLENNIUM® UNIT AND 1363-DD.

1367 OERTLI® VFI ADAPTOR

FOR CONNECTION BETWEEN OERTLI® UNIT AND 1363-DD.

1368 GEUDER® VFI ADAPTOR

FOR CONNECTION BETWEEN GEUDER® UNIT AND 1363-DD.

1369 RUCK® VFI ADAPTOR

FOR CONNECTION BETWEEN RUCK® UNIT AND 1363-DD.

NOTE: 1363-DD FITS ACCURUS® UNIT WITHOUT ADAPTOR


ACCESSORIES FOR VFE

1362-D2 DISPOSABLE VISCOUS FLUID EXTRACTION KIT

- VFE KIT IN COMBINATION WITH TWO 10CC SYRINGES FOR BIMANUAL ASPIRATION
- REDUCES ASPIRATION TIME BY 50%
- STERILE, BOX/5

PATIENT RECORDS LABEL

PATIENT RECORDS LABEL FOR OPTIMAL TRACEABILITY FOUR LABELS ARE SUPPLIED WITH EACH VARIETY.

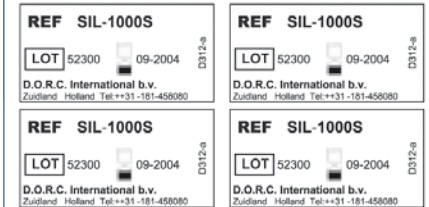

For additional information please contact:

D.O.R.C. International B.V.
Scheijdelweg 2, 3214 VN Zuidland
The Netherlands
Phone: +31 181 45 80 80
Fax: +31 181 45 80 90
e-mail: sales@dorc.nl

In Latin America please contact:

D.O.R.C. Latin America
Scheijdelweg 2, 3214 VN Zuidland
The Netherlands
Phone: +1 305-663-5662
Fax: +1 305-663-7343
e-mail: la-sales@dorc.nl

In Middle East please contact:

D.O.R.C. Middle East, Africa and India
PO BOX 212218, Dubai, UAE
Phone: ++971 (0) 50 911 65 87
Fax: ++31 (0) 181 45 80 90
e-mail: cvds@dorc.nl

In Scandinavia please contact:

D.O.R.C. Scandinavia AB
Bettorpsgatan 10, SE-701 11 Örebro, Sweden
Phone: +46 19121514
Fax: +46 19188204
e-mail: order@dorc.se

In France please contact:

D.O.R.C. France S.A.R.L.
30, Boulevard Pasteur, 75015 Paris, France
Phone: +33 1 53 58 14 73
Fax: +33 1 53 58 14 75
e-mail: ventes@dorc.fr

In Germany please contact:

D.O.R.C. Deutschland GmbH
Charlottenstraße 68, 10117 Berlin, Deutschland
Phone: +49 30 201 88364
Fax: +49 30 201 88365
e-mail: bestellung@dorc.de

AUTHORIZED D.O.R.C. DEALER

