



FOCUS ON PRECISION

FOCUS ON INNOVATION

FOCUS ON QUALITY



FOCUS ON
XENON BRIGHTSTAR

Xenon BrightStar ILLUMINATION SYSTEM

THE D.O.R.C. XENON BRIGHT STAR ILLUMINATION SYSTEM PROVIDES THE NEW GENERATION ULTIMATE BRIGHT WHITE ILLUMINATION FOR IMPROVED VISIBILITY DURING SURGERY, EVEN WITH THE SMALLEST LIGHT FIBERS!

D.O.R.C. HAS ALSO RESPONDED TO CLINICAL CONCERNS REGARDING MINIMIZING THE INCIDENCE OF PHOTO-TOXICITY, ESPECIALLY WHEN SURGICAL CONDITIONS DICTATE WORKING CLOSE TO THE RETINA.

THE D.O.R.C. XENON BRIGHT STAR IS THE FIRST ILLUMINATION SYSTEM TO OFFER A SELECTION OF ULTRAVIOLET CUT-OFF FILTERS, GIVING THE SURGEON ULTIMATE SAFETY CONTROL OVER THE SPECTRUM OF LIGHT RAYS EXPOSED TO THE RETINA.

COUPLED WITH AN EXTENSIVE ARRAY OF D.O.R.C. ILLUMINATION ACCESSORIES, XENON BRIGHT STAR PROVIDES THE ULTIMATE IN PERFORMANCE, VERSATILITY AND VALUE.



SAFETY:

- FEATURING MULTIPLE SAFETY FILTER SELECTION
- HIGH INTENSIVE ILLUMINATION IN COMBINATION WITH INCREASED SAFETY

SAFETY

D.O.R.C. IS THE FIRST MANUFACTURER TO INCLUDE A SELECTION OF ULTRA-VIOLET CUT-OFF FILTERS WHICH GIVES THE SURGEON PRECISE CONTROL OF THE ILLUMINATION SPECTRUM NEEDED FOR VARIOUS SURGICAL SITUATIONS. THE FOUR FILTER CUT-OFF CHOICES INCLUDE:

- 420NM CUT-OFF: STANDARD
- 435NM CUT-OFF: FOR CORE VITRECTOMY AND GENERAL MEMBRANE REMOVAL
- 475NM CUT-OFF: FOR REMOVAL OF MEMBRANES ADHERENT TO THE RETINA
- 515NM CUT-OFF: FOR MACULAR REPAIR

THESE FILTERS PROTECT THE RETINA AND ALLOWS WORKING CLOSE TO THE RETINA FOR LONGER TIMES.

HIGHLIGHTS

- HIGH INTENSIVE OUTPUT IN COMBINATION WITH SAFE / MINIMAL TOXIC LIGHT
- LIGHT TIME INDICATOR
- HIGH INTENSIVE XENON BULB VERTICALLY POSITIONED WITH A LONG LIFE TIME DUTY CYCLE
- EASY EXCHANGEABLE XENON BULB
- FOUR SELECTABLE UV CUT-OFF FILTERS!

For additional information please contact:

D.O.R.C. International B.V.
Scheijdelweg 2, 3214 VN Zuidland
The Netherlands
Phone: +31 181 45 80 80
Fax: +31 181 45 80 90
E-mail: sales@dorc.nl

In US / Canada please contact:

Dutch Ophthalmic USA
10 Continental Drive
Bldg 1, Exeter, NH 03833, U.S.A.
Phone: +1 800-75-DUTCH or
+1 603-778-6929
Fax: +1 603-778-0911
E-mail: sales@dutchophthalmicusa.com

In Latin America please contact:

D.O.R.C. Latin America
Scheijdelweg 2, 3214 VN Zuidland
The Netherlands
Phone: +1 305-663-5662
Fax: +1 305-663-7343
E-mail: la-sales@dorc.nl

In Middle East please contact:

D.O.R.C. Middle East, Africa and India
PO BOX 212218, Dubai, United Arab Emirates
Phone: 00971 48857613
Fax: 00971 48857614
E-mail: cvds@dorc.nl

In Scandinavia please contact:

D.O.R.C. Scandinavia AB
Bettorpsgatan 10, SE-701 11 Örebro, Sweden
Phone: +46 19121514
Fax: +46 19188204
E-mail: order@dorc.se


In France please contact:

D.O.R.C. France S.A.R.L.
30, Boulevard Pasteur, 75015 Paris, France
Phone: +33 1 53 58 14 73
Fax: +33 1 53 58 14 75
E-mail: ventes@dorc.fr

In Germany please contact:

D.O.R.C. Deutschland GmbH
Charlottenstraße 68, 10117 Berlin, Deutschland
Phone: +49 30 201 88364
Fax: +49 30 201 88365
E-mail: bestellung@dorc.de

AUTHORIZED D.O.R.C. DEALER

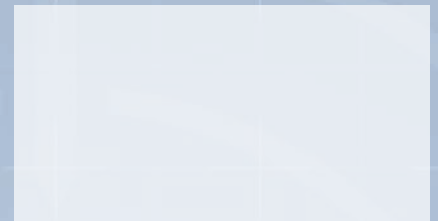

FOCUS MORE AT WWW.DORC.NL